




TYCOON SWISS MEDICAL

Most technologically advanced
device for Cryo neuromodulation

CRYO-S PAINLESS

Recommended for Pain Management specialists:


- Anesthesiologists
- Neurosurgeons
- Orthopedists
- Spine surgeons
- Interventional Radiologist

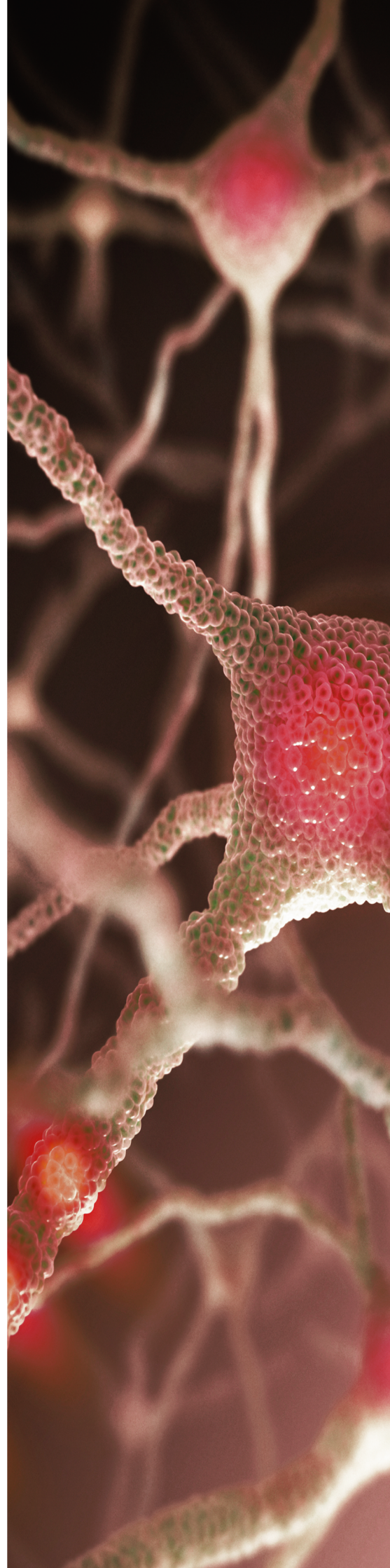


Peripheral nerve structure




1. Peripheral neuropathy

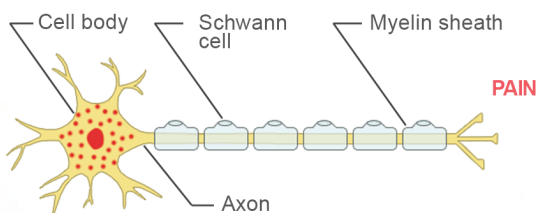


2. Interventional cryo neuromodulation

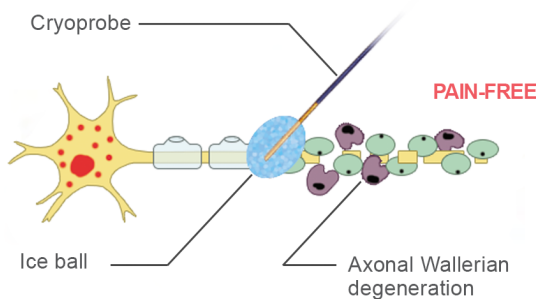

3. Axonal nerve recovery process

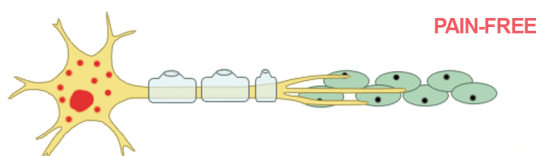

4. Complete nerve recovery

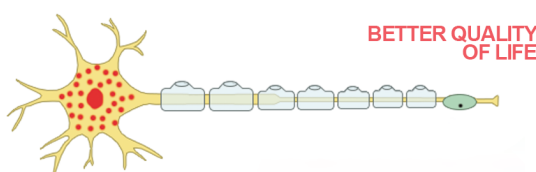

Life without pain

Cryo neuromodulation is a therapeutic method based on the temporary interruption of sensory functions in selected structures of the nervous system by application of low temperature. Cryo neuromodulation is gaining more acceptance as an innovative method in pain relief. It uses the process of analgesia, during which the ice crystals created by the cryosurgical system destroy the elements of the nerve tissue carrying pain information. Cryo neuromodulation does not damage nerve structures permanently which is why nerve tissues can regenerate slowly with no risk of post-procedural neuroma. Cryo neuromodulation is a minimally invasive and safe procedure recommended especially where traditional methods prove to be unsatisfactory.

Mechanism of cryo neuromodulation

Large myelinated sensory or mixed nerves responsible for peripheral pain. Best patients for cryo neuromodulation are those with chronic pain at the level of over 5 in visual analogue scale (VAS).

Axonotmesis - disruption of nerve cell axon with Wallerian degeneration occurring below and slightly proximal to the site of the injury. Axons and their myelin sheath are damaged but Schwann cells, endoneurium, perineurium and epineurium remain intact.

The rate of outgrowth of regenerating nerve fibre is about 1 mm per day, so that the final recovery of conduction to the target structure depends on the quality of freezing. Even if the peripheral nerve structures recover after some weeks, the functionality of the sensory nerve still remains blocked. In mixed nerves, mobility functions recover much faster than in sensory. Cryo neuromodulation is also a recommended and safe method for pain of mixed nerves' origin (e.g. "frozen shoulder").

Complete sensory nerve recovery, including its functionality, takes place within 6-12 months after cryo neuromodulation improving patient's quality of life - life without pain.

Clinical application

- Craniofacial pain secondary to trigeminal neuralgia, posterior auricular neuralgia, and glossopharyngeal neuralgia.
- Chest wall pain with multiple conditions including post-thoracotomy neuromas, persistent pain after rib fractures and post herpetic neuralgia in thoracic distribution.
- Abdominal and pelvic pain secondary to ilioinguinal, iliohypogastric, genitofemoral, subgastric neuralgia, pudendal neuralgia.
- Low back pain and lower extremity pain secondary to lumbar facet joint pathology, pseudosciatica.
- Pain involving intraspinous ligament of supragluteal nerve.
- Sacroiliac joint pain, cluneal neuralgia, obturator neuritis.
- Knee pain prior or secondary to infrapatellar branch of the saphenous nerve, anterior femoral cutaneous nerve, genicular nerve after total knee arthroplasty.
- Various types of peripheral neuropathy.
- Trigger points.
- Phantom pain.
- Upper extremity pain secondary to suprascapular neuritis and other conditions of peripheral neuritis.
- Diabetic neuropathy.
- Morton's neuroma, heel pain.

The advantages of cryo neuromodulation

- No neuroma formation - no risk of secondary pain.
- High efficiency: pain reduction from 6 months to 2 years.
- No scar tissue formation.
- Suitable for patients with pacemakers and stimulators.
- No risk of vessel proliferation and obliteration.
- Can be repeated - nerve grows back.
- Simple diagnostics: fluoroscopy or echography.
- Percutaneous procedure, microinvasive under local anesthesia.
- Fast return to normal activity - no hospitalization.
- Can be performed under USG or X-ray.




Great occipital nerve
(headic pain)



Suprascapular nerve
(frozen shoulder)



Genicular nerve
(knee pain)

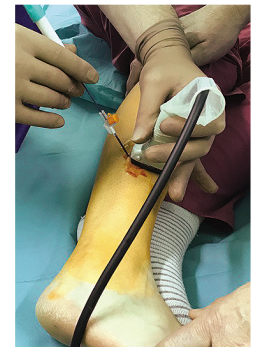



Painful neuroma



Morton's neuroma (foot pain)



Genitofemoral nerve (pelvic pain)

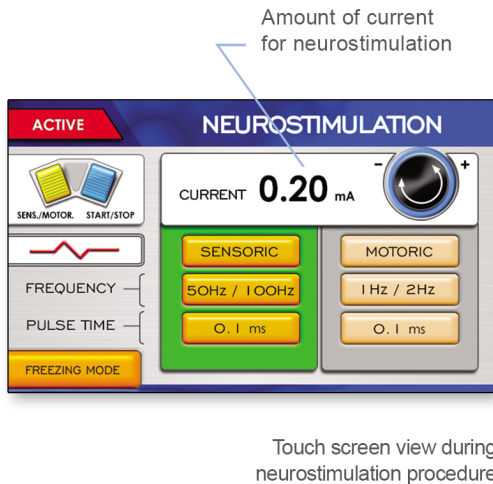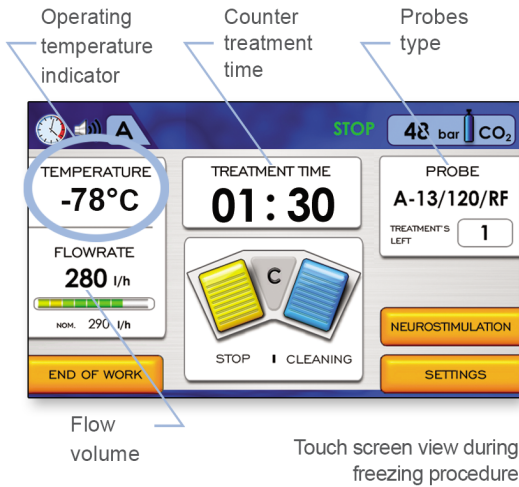

Grounding pad

Cooling source

CRYO-S PAINLESS state-of-the-art cryosurgical device manufactured by METRUM CRYOFLEX is the next generation of apparatus used by many experts in the field since 1992.

The working medium for CRYO-S PAINLESS is carbon dioxide: CO₂ (-78°C), very efficient and easy to use gas. A 10 liter cylinder of CO₂ serves for approximately 60 procedures.

CRYO-S PAINLESS is controlled by a microprocessor and all the parameters are displayed and monitored on LCD screen.

Neurostimulation

Diagnostic nerve stimulation built in CRYO-S PAINLESS device helps to distinguish between sensory and motor nerves. It prevents from accidental freezing of motor nerves, which is a reversible effect, however the possible discomfort connected with temporary paresis may entirely challenge pain therapy.

Integrated neurostimulator provides two types of stimulation (sensory and motor) and a variety of options for manual settings (without freezing).

Stimulation frequency setting:

- For motor stimulation: 1 and 2 Hz.
- For sensory stimulation: 50, 100, 150, 200 Hz.
- Impulse setting: 0.1, 0.2, 0.5, 1.0, 2.0 – 5.0 ms.

Grounding pad for use with CRYO-S PAINLESS device

Single use grounding pad dedicated for CRYO-S PAINLESS is easy to use and provides more precise neurostimulation than any other grounding pad. Available accessories including connector cables, grounding pads are always available with CRYO-S PAINLESS.

CRYO-S PAINLESS device

- Mode selection probe, cleaning and freezing can be performed automatically using footswitch or touch screen which allows to keep the site of a procedure under sterile conditions.
- Electronic communication (chip system) between the connected probe and device allows recognition of optimal operating parameters and auto-configures to cryoprobe characteristics. Pressure and gas flow are set automatically, any manual adjustments not necessary.
- Contains automatic cryoprobe test system and two freezing modes (continuous and sequential). Probe test system includes autocleaning and short freezing test.
- Cryoprobe temperature, cylinder pressure, gas flow inside of cryoprobe and procedure time are displayed during freezing test.
- Built-in voice communication.
- Built-in neurostimulation (sensory, motor).

Operated with footswitch
or touch screen



TOUCH

Electronic communication
between the probe and device



RFID

Two freezing modes



MODE

Voice communication



VOICE


Neurostimulation
(sensory, motor)





PNS



CE 2274

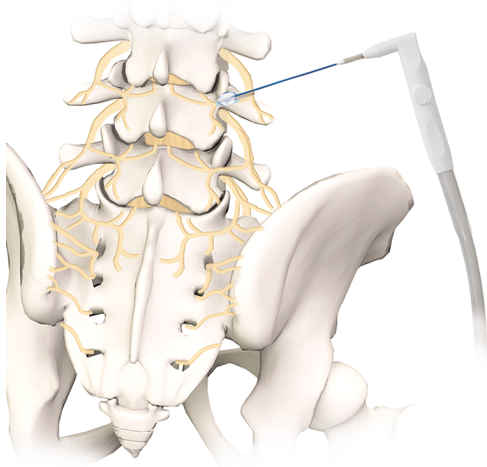

CE 2274

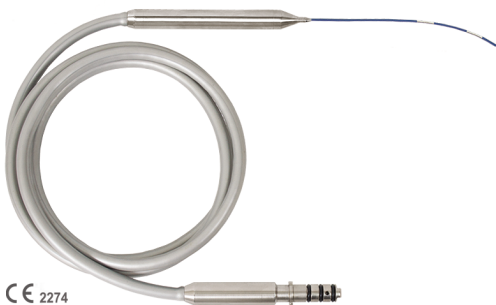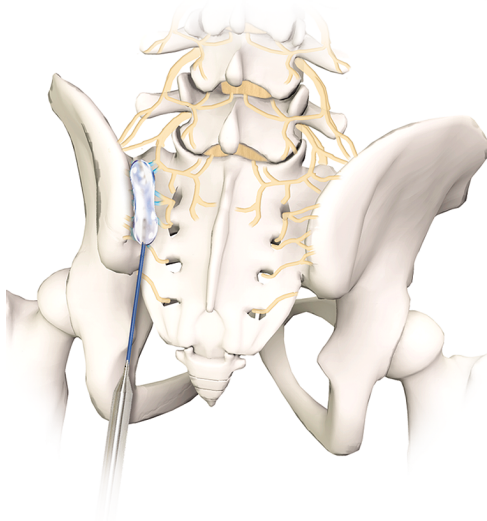

CE 2274

Angular cryoprobe

The largest selection of cryoprobes - the thinnest cryoprobes for cryo neuromodulation available on the global pain management market.

We offer both single use and reusable cryoprobes which are easy to sterilize (steam, gas). The special construction enables precise freezing and fast defrosting. The tip of a probe inserted percutaneously in sharp or blunt. Cryoprobes that aim to stick to the nerve more closely are rounded and inserted through disposable or reusable cannulas (e.g. angiocatheters).

Apart from the tip, cryoprobes are covered with medical teflon which allows for stimulation at the probe tip only.

- Gauge from 21 up to 14.
- Angular, autoclavable cryoprobes.
- Convenient handle of the cryoprobe ensures safety of the procedures conducted with fluoroscopic guidance
- More safety - the operator's hand is not exposed to the harmful X-ray radiation.


Sacroiliacjoint cryoprobe (SIJ)

New SIJ (sacroiliacjoint) cryoprobe represents a new vision of performing cryo neuromodulation for chronic lower back pain, including ablation of the sacroiliac region. With independent active areas in a single probe, procedure time is generally reduced.

- Design is intended to create a continuous cryoablation with no gaps (similar to RF palisade method).
- A single insertion point reduces time of procedure.
- Designed to lesion the peripheral nerves at S1-S3.
- Ergonomic handle enhances maneuverability.
- Introducer not needed.

Single-use cryoprobe

recommended for USG, CT, X-ray scan

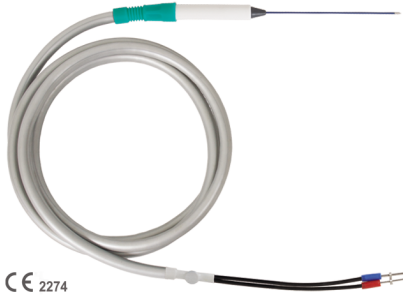

CE 2274

| | | | | | | |
|-----------------------|-------------------|------------|-------------------------|------------|-------------------|----------------|
| CO₂ | STERILE EO | 2 | DOUBLE PACKAGING | PNS | MICRO CHIP | 4 YEARS |
| GAS TYPE | READY TO USE | SINGLE USE | DOUBLE PACKAGING | PNS | MICRO CHIP | STORAGE TIME |

REF. SN06080921 Needle tip***, \varnothing 0.8 mm (21Ga), Length 90mm


REF. SN06131200 Triangular tip*, \varnothing 1.3 mm (18Ga), Length 100mm

REF. SN06300200 Triangular tip*, \varnothing 2.0 mm (14Ga), Length 120mm

REF. SJA Multi-use adapter for single-use cryoprobe

Reusable straight cryoprobe

recommended for USG, or X-ray scan

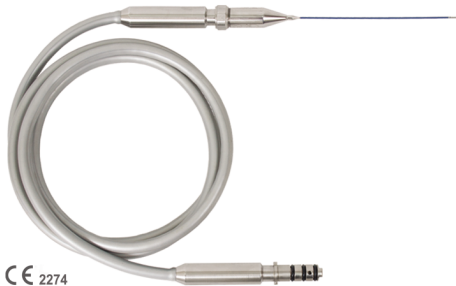

CE 2274

| | | | | | | |
|-----------------------|----------------------|----------------------|--------------------------|------------|-------------------|----------------|
| CO₂ | AUTO CLAVE | ETO GAS | 50 | PNS | MICRO CHIP | 3 YEARS |
| GAS TYPE | STERILISATION METHOD | STERILISATION METHOD | NUMBER OF STERILISATIONS | PNS | MICRO CHIP | WORKING LIFE |

REF. N06080921 Needle tip***, \varnothing 0.8 mm (21Ga), Length 90 mm

REF. N06080421 Needle tip***, \varnothing 0.8 mm (21Ga), Length 40 mm

REF. N06131201 Triangular tip*, \varnothing 1.3 mm (18Ga), Length 100 mm

REF. N06131211 Round tip**, \varnothing 1.3 mm (18Ga), Length 100 mm

REF. N06300201 Triangular tip*, \varnothing 2.0 mm (14Ga), Length 120mm

Reusable angular cryoprobe

recommended for X-ray or CT scan



CE 2274

| | | | | | | |
|-----------------------|----------------------|----------------------|--------------------------|------------|-------------------|----------------|
| CO₂ | AUTO CLAVE | ETO GAS | 50 | PNS | MICRO CHIP | 3 YEARS |
| GAS TYPE | STERILISATION METHOD | STERILISATION METHOD | NUMBER OF STERILISATIONS | PNS | MICRO CHIP | WORKING LIFE |


REF. N16300200 Triangular tip*, \varnothing 2.0 mm (14Ga), Length 120mm

REF. N16131200 Triangular tip*, \varnothing 1.3 mm (18Ga), Length 100mm

REF. N16131210 Round tip**, \varnothing 1.3 mm (18Ga), Length 100 mm

Reusable SIJ cryoprobe

recommended for X-ray, USG or CT scan

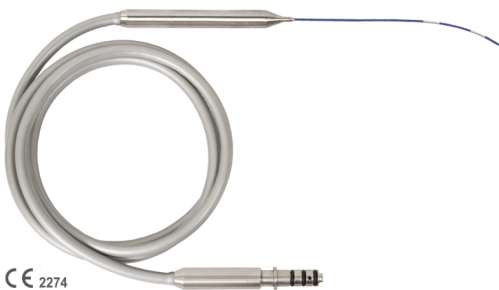

CE 2274

| | | | | | | |
|-----------------------|----------------------|----------------------|--------------------------|------------|-------------------|----------------|
| CO₂ | AUTO CLAVE | ETO GAS | 25 | PNS | MICRO CHIP | 3 YEARS |
| GAS TYPE | STERILISATION METHOD | STERILISATION METHOD | NUMBER OF STERILISATIONS | PNS | MICRO CHIP | WORKING LIFE |

REF. N06302200 Triangular tip*, \varnothing 2.0 mm (14Ga), Length 120 mm

| | | |
|---------------------|-----------------|-------------------|
| * Triangular tip | ** Round tip | *** Needle tip |
|---------------------|-----------------|-------------------|



TYCOON SWISS MEDICAL SA
Via San Salvatore n°10 CH-6902 Paradiso
Telephone: +41 91 993 0504
E-mail: info@tycoonsm.com

MANUFACTURER

M C METRUM CRYOFLEX

METRUM CRYOFLEX Sp. z o.o., Sp.K.

Production
ul. Zielna 29
05-082 Blizne Łaszczyńskiego
POLAND, EU

Headoffice
ul. Kolejowa 16A
05-092 Łomianki
POLAND, EU

Tel: +4822 33 13 750 office@metrum.com.pl
+4822 33 13 830 <http://www.metrum.com.pl>
Fax: +4822 33 13 766